


BASIC MASS SPEC BENCH 72W

WITH DUAL PUMP ENCLOSURE



OL-MS-BAS72-D

Overall Bench Dimensions: 72W x 33D x 37H

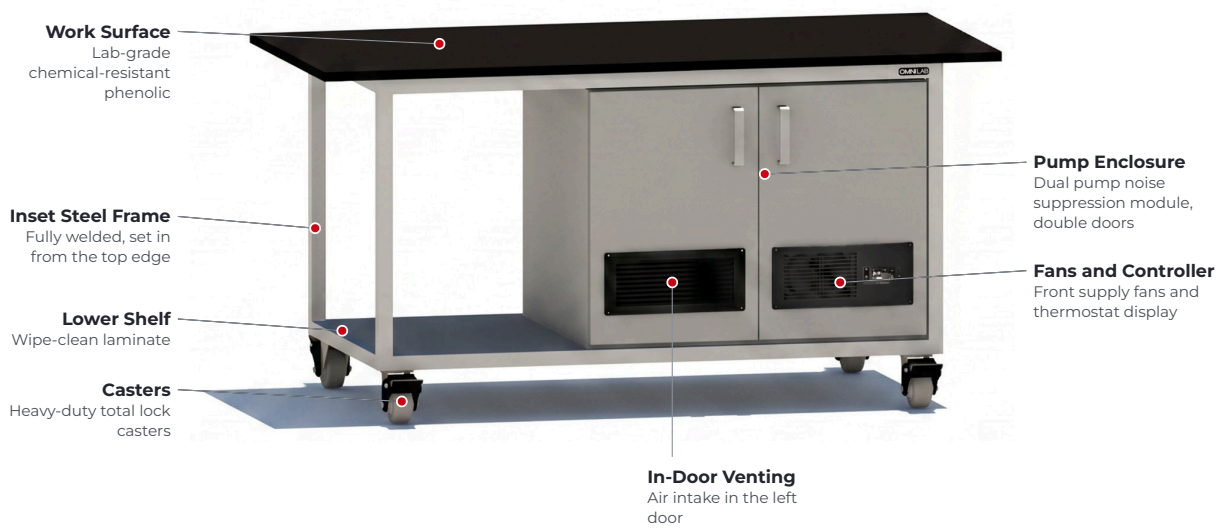

Open Mass Spec Bench with lab-grade chemical-resistant phenolic work surface and inset steel frame with lower shelf on heavy-duty total lock casters. Includes a noise suppression module for dual vacuum pumps with acoustic foam paneling, push-pull cooling (440 CFM combined rated exhaust), in-door venting, double doors, five rear hose grommets, and two oil drip pans with shock absorbers. Also includes one 120V 8-outlet surge suppressor with a 12-foot cord. Overall bench size: 72W x 33D x 37H (Pump Enclosure Interior: 32-1/2W x 25-3/4D x 24H).

SPECIFICATIONS

OVERALL SIZE 72W x 33D x 37H	WORK SURFACE 1" lab-grade chemical-resistant phenolic	FRAME Inset fully welded steel, epoxy powder coat finish, with wipe-clean laminate shelf
PUMP ENCLOSURE Acoustic foam lined, 440 CFM push-pull cooling, in-door venting, thermostat controller	LOAD CAPACITY 1,200 lb uniform load	CASTERS Heavy-duty total lock

BASIC MASS SPEC BENCH 72W

WITH DUAL PUMP ENCLOSURE


OMNILAB

SOLUTIONS

COMPONENTS INCLUDED

QTY	MODEL #	DESCRIPTION
1	OL3IT72	Open Instrument Bench 72W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Inset Steel Frame (60W x 31D x 30-1/4H) with Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity
1	OL3NS-D	Bench-Mounted Noise Suppression Module for Dual Vacuum Pumps 36W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with In-Door Venting, Front Supply Fans and Quad Rear Exhaust Fans (440 CFM Combined Rated Exhaust), Five (5) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Double Doors
2	OL3ODP-XL	Oil Drip Pan for Pump with Mounted Shock Absorbers, Large
1	OL3E	Tripp Lite® Isobar® Surge Suppressor, Eight (8) 120V NEMA 5-15R Outlets in Isolated Filter Banks for EMI/RFI Noise Protection, 3,840 Joule Rating, 12 Amp Resettable Circuit Breaker (1,440 Watt Capacity), UL Listed, Diagnostic LEDs, All-Metal Housing, 12 Ft Cord, Bench Mounted

PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

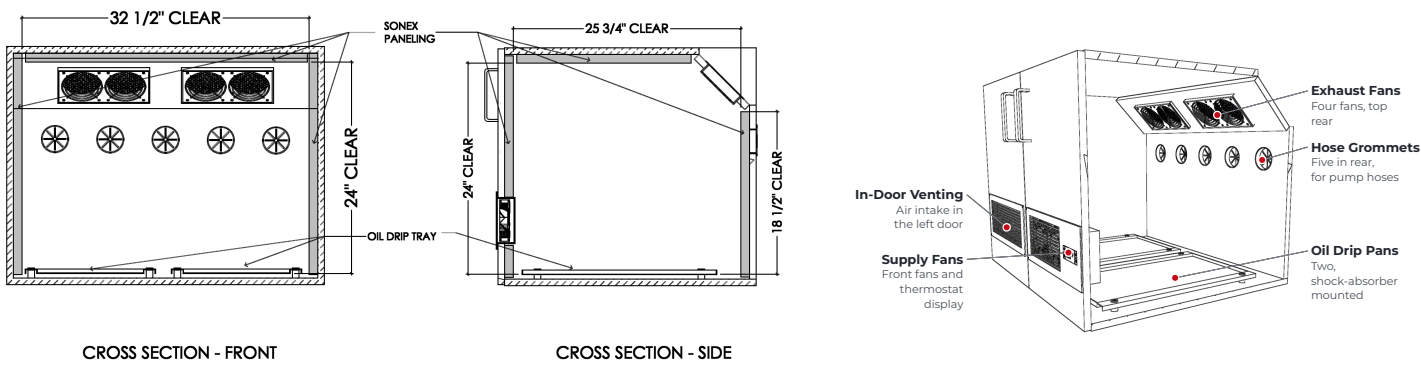

Interior clear dimensions: 32-1/2W x 25-3/4D x 24H (18-1/2H at the rear wall, below the angled exhaust fans). Confirm pumps, hoses, cords and service clearance fit before ordering.

FINISH OPTIONS



Main image shown in **Gray** (standard). **Also available in White** (shown here). Additional frame colors may be available upon request.