

TOTAL LCMS BENCH

WITH SINGLE PUMP ENCLOSURE



OL-MS-LCMS-TOT-S

Overall Bench Dimensions: 72W x 33D x 37H



Enclosed LCMS Bench with lab-grade chemical-resistant phenolic work surface and inset steel frame with lower shelf on heavy-duty total lock casters. Includes a noise suppression module for a single vacuum pump with acoustic foam paneling, push-pull cooling (220 CFM combined rated exhaust), a vented door, three rear hose grommets, a slide-out maintenance shelf with 30 degree tilt for oil drainage, and an oil drip pan with shock absorbers; an enclosed waste container cabinet with a vision panel in the door, top-front waste line access, and a slide-out shelf; and a drawer bank with three full-extension drawers and a pull-out writing shelf. Also includes a sit-stand PC workstation (monitor arm and keyboard tray) with a vertical CPU/UPS holder, plus two 120V 8-outlet surge suppressors with 12-foot cords. Overall bench size: 72W x 33D x 37H (Pump Enclosure Interior: 18-1/2W x 25-3/4D x 22-1/4H).

SPECIFICATIONS

OVERALL SIZE 72W x 33D x 37H	WORK SURFACE 1" lab-grade chemical-resistant phenolic	FRAME Inset fully welded steel, epoxy powder coat finish, with wipe-clean laminate shelf
PUMP ENCLOSURE Acoustic foam lined, 220 CFM push-pull cooling, thermostat controller	CASTERS Heavy-duty total lock	STORAGE AND PANELS Wipe-clean laminate in a color that coordinates with the frame

TOTAL LCMS BENCH



WITH SINGLE PUMP ENCLOSURE



COMPONENTS INCLUDED

QTY	MODEL #	DESCRIPTION
1	OL3IT72E	Enclosed Instrument Bench 72W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Inset Steel Frame (60W x 31D x 30-1/4H) with Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity
1	OL3NS-S-R	Bench-Mounted Noise Suppression Module for Single Vacuum Pump 22W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with Dual Front Supply Fans and Dual Rear Exhaust Fans (220 CFM Combined Rated Exhaust), Three (3) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Door Hinged Right
1	OL3SOS-S	Slide-Out Pump Maintenance Shelf with 30 Degree Tilt for Above Noise Suppression Module, Allows Gravity-Assisted Oil Drainage for Pump Oil Changes, Single Pump Shelf
1	OL3ODP-XL	Oil Drip Pan for Pump with Mounted Shock Absorbers, Large
1	OL3DB16	Drawer Bank 15-3/4W x 22D x 27H with Three (3) Full Extension Lab Drawers and One (1) Pull-Out Utility Shelf
1	OL3WS21	Enclosed Waste Container Storage Cabinet 20-3/4W x 28-1/2D x 27H with Vision Panel in Door and Waste Line Access from Top Front, Door Hinged Left
1	OL3WSOS	Open Slide-Out Shelf for Waste Container Storage with Chemical Resistant Laminate Surface 16W x 22D x 1-1/2H
1	OL3PCW-C	Ergotron® Work Surface Mounted Sit-Stand Workstation, Ergonomic Monitor Arm and Keyboard Tray with Independent Height Adjustment, 20 Inch Height Range, Holds Monitors up to 42 Inch, Slide-Out Left or Right Mouse Tray, Wrist Rest, Integrated Cable Management, Clamp or Grommet Mount, Clamps to Surfaces up to 2.5 Inch Thick (Including Work Surface and Frame)
1	OL3CPU	Ergotron® Side-Mounted Vertical CPU/UPS Holder, Adjustable to Fit Units 2 to 8 Inch Wide, 40 lb (18 kg) Capacity, Open Design for Access to Ports and Power Buttons, Steel and Aluminum Construction
2	OL3E	Tripp Lite® Isobar® Surge Suppressor, Eight (8) 120V NEMA 5-15R Outlets in Isolated Filter Banks for EMI/RFI Noise Protection, 3,840 Joule Rating, 12 Amp Resettable Circuit Breaker (1,440 Watt Capacity), UL Listed, Diagnostic LEDs, All-Metal Housing, 12 Ft Cord, Bench Mounted

PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

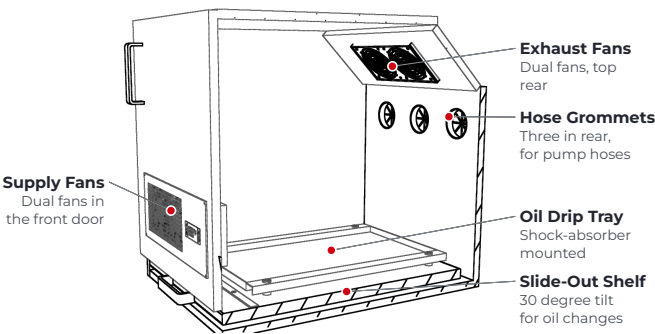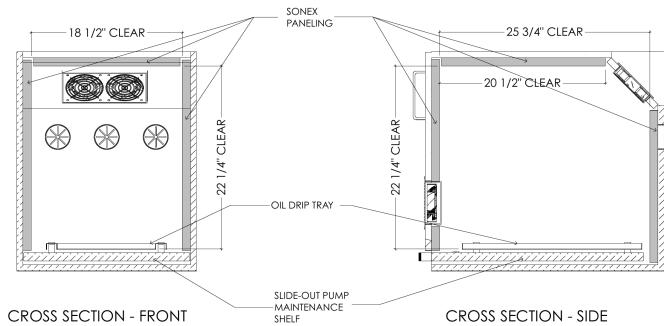

Interior clear dimensions: 18-1/2W x 25-3/4D x 22-1/4H (20-1/2D at the top, where the exhaust fans sit at an angle). Confirm pump, hoses, cords and service clearance fit before ordering.

FINISH OPTIONS



Main image shown in **Gray** (standard). Also available in **White** (shown here). Additional frame colors may be available upon request.