

# CORE GCMS BENCH

WITH DUAL PUMP ENCLOSURE



OL-MS-GCMS-CORE-D

Overall Bench Dimensions: 72W x 33D x 37H



Enclosed GCMS Bench with lab-grade chemical-resistant phenolic work surface and inset steel frame with lower shelf on heavy-duty total lock casters. Includes a noise suppression module for dual vacuum pumps with acoustic foam paneling, push-pull cooling (440 CFM combined rated exhaust), in-door venting, double doors, five rear hose grommets, and two oil drip pans with shock absorbers; and a drawer bank with three full-extension drawers and a pull-out writing shelf. Also includes two 120V 8-outlet surge suppressors with 12-foot cords. Overall bench size: 72W x 33D x 37H (Pump Enclosure Interior: 32-1/2W x 25-3/4D x 24H).

SPECIFICATIONS

|                                                                                                          |                                                          |                                                                                             |
|----------------------------------------------------------------------------------------------------------|----------------------------------------------------------|---------------------------------------------------------------------------------------------|
| OVERALL SIZE<br>72W x 33D x 37H                                                                          | WORK SURFACE<br>1" lab-grade chemical-resistant phenolic | FRAME<br>Inset fully welded steel, epoxy powder coat finish, with wipe-clean laminate shelf |
| PUMP ENCLOSURE<br>Acoustic foam lined, 440 CFM push-pull cooling, in-door venting, thermostat controller | CASTERS<br>Heavy-duty total lock                         | STORAGE AND PANELS<br>Wipe-clean laminate in a color that coordinates with the frame        |

# CORE GCMS BENCH


WITH DUAL PUMP ENCLOSURE



## COMPONENTS INCLUDED

| QTY | MODEL #   | DESCRIPTION                                                                                                                                                                                                                                                                                                                                                                                                                                         |
|-----|-----------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1   | OL3IT72E  | Enclosed Instrument Bench 72W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Inset Steel Frame (60W x 31D x 30-1/4H) with Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity                                                                                                                                                                                                                     |
| 1   | OL3DB21   | Drawer Bank 20-3/4W x 22D x 27H with Three (3) Full Extension Lab Drawers and One (1) Pull-Out Utility Shelf                                                                                                                                                                                                                                                                                                                                        |
| 1   | OL3NS-D   | Bench-Mounted Noise Suppression Module for Dual Vacuum Pumps 36W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with In-Door Venting, Front Supply Fans and Quad Rear Exhaust Fans (440 CFM Combined Rated Exhaust), Five (5) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Double Doors |
| 2   | OL3ODP-XL | Oil Drip Pan for Pump with Mounted Shock Absorbers, Large                                                                                                                                                                                                                                                                                                                                                                                           |
| 2   | OL3E      | Tripp Lite® Isobar® Surge Suppressor, Eight (8) 120V NEMA 5-15R Outlets in Isolated Filter Banks for EMI/RFI Noise Protection, 3,840 Joule Rating, 12 Amp Resettable Circuit Breaker (1,440 Watt Capacity), UL Listed, Diagnostic LEDs, All-Metal Housing, 12 Ft Cord, Bench Mounted                                                                                                                                                                |

## PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

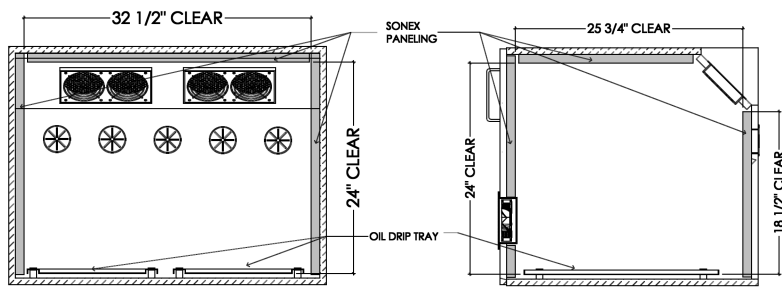

CROSS SECTION - FRONT

CROSS SECTION - SIDE

**Interior clear dimensions: 32-1/2W x 25-3/4D x 24H** (18-1/2H at the rear wall, below the angled exhaust fans). Confirm pumps, hoses, cords and service clearance fit before ordering.

## FINISH OPTIONS



Main image shown in **Gray** (standard). **Also available in White** (shown here). Additional frame colors may be available upon request.

