

ULTIMATE GCMS BENCH

WITH SINGLE PUMP ENCLOSURE



OL-MS-GCMS-ULT-S

Overall Bench Dimensions: 96W x 33D x 37H



Two-Piece Enclosed GCMS Bench (60W MS Bench and 36W PC/Ancillary Instrument Bench) with lab-grade chemical-resistant phenolic work surfaces and inset steel frames with lower shelves on heavy-duty total lock casters. Includes a noise suppression module for a single vacuum pump with acoustic foam paneling, push-pull cooling (220 CFM combined rated exhaust), a vented door, three rear hose grommets, a slide-out maintenance shelf with 30 degree tilt for oil drainage, and an oil drip pan with shock absorbers; and a drawer bank with three full-extension drawers and a pull-out writing shelf. Also includes a PC accessory kit with a monitor arm, a sit-stand keyboard arm, and a vertical CPU/UPS holder, plus two 120V 8-outlet surge suppressors with 12-foot cords. Overall bench size: 96W x 33D x 37H (Pump Enclosure Interior: 18-1/2W x 25-3/4D x 22-1/4H).

SPECIFICATIONS

OVERALL SIZE 96W x 33D x 37H	WORK SURFACE 1" lab-grade chemical-resistant phenolic	FRAME Fully welded steel, epoxy powder coat finish, with wipe-clean laminate shelves
PUMP ENCLOSURE Acoustic foam lined, 220 CFM push-pull cooling, thermostat controller	CASTERS Heavy-duty total lock	STORAGE AND PANELS Wipe-clean laminate in a color that coordinates with the frame

ULTIMATE GCMS BENCH



WITH SINGLE PUMP ENCLOSURE



COMPONENTS INCLUDED

QTY	MODEL #	DESCRIPTION
1	OL3IT60E	Enclosed Instrument Bench 60W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Inset Steel Frame (48W x 31D x 30-1/4H) with Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity
1	OL3IT36E	Enclosed Instrument Bench 36W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Inset Steel Frame (30W x 31D x 30-1/4H) with Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity
1	OL3DB21	Drawer Bank 20-3/4W x 22D x 27H with Three (3) Full Extension Lab Drawers and One (1) Pull-Out Utility Shelf
1	OL3NS-S-R	Bench-Mounted Noise Suppression Module for Single Vacuum Pump 22W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with Dual Front Supply Fans and Dual Rear Exhaust Fans (220 CFM Combined Rated Exhaust), Three (3) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Door Hinged Right
1	OL3SOS-S	Slide-Out Pump Maintenance Shelf with 30 Degree Tilt for Above Noise Suppression Module, Allows Gravity-Assisted Oil Drainage for Pump Oil Changes, Single Pump Shelf
1	OL3ODP-XL	Oil Drip Pan for Pump with Mounted Shock Absorbers, Large
1	OL3MA-C	Ergotron® Work Surface Mounted Ergonomic Monitor Arm, Clamp Mount, Holds Monitors up to 34 Inch and 22 lb, VESA 75 x 75 and 100 x 100 Compatible, Height, Tilt, Swivel and Rotation Adjustment, Integrated Cable Management
1	OL3AKMP-SS	Ergotron® Frame Mounted Sit-Stand Ergonomic Keyboard Arm, 13 Inch Height Adjustment, Back-Tilt Keyboard Tray for Neutral Wrist Position, Slide-Out Mouse Tray, Full Swivel and Extension, Folds Upright for Storage, Built-In Rotation Stop
1	OL3CPU	Ergotron® Side-Mounted Vertical CPU/UPS Holder, Adjustable to Fit Units 2 to 8 Inch Wide, 40 lb (18 kg) Capacity, Open Design for Access to Ports and Power Buttons, Steel and Aluminum Construction
2	OL3E	Tripp Lite® Isobar® Surge Suppressor, Eight (8) 120V NEMA 5-15R Outlets in Isolated Filter Banks for EMI/RFI Noise Protection, 3,840 Joule Rating, 12 Amp Resettable Circuit Breaker (1,440 Watt Capacity), UL Listed, Diagnostic LEDs, All-Metal Housing, 12 Ft Cord, Bench Mounted

PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

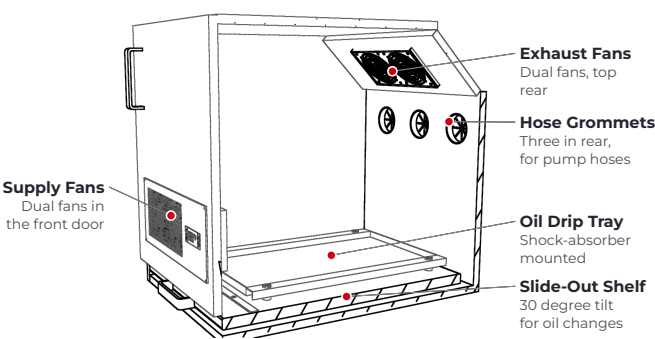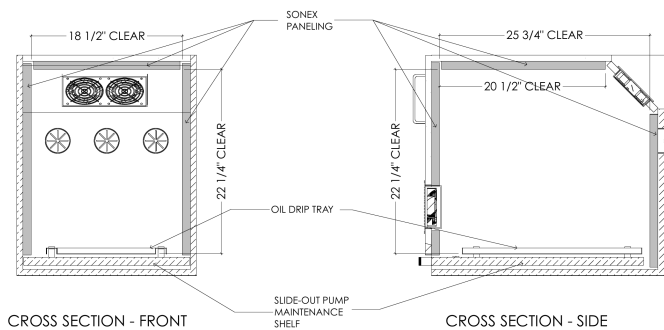

Interior clear dimensions: 18-1/2W x 25-3/4D x 22-1/4H (20-1/2D at the top, where the exhaust fans sit at an angle). Confirm pump, hoses, cords and service clearance fit before ordering.

FINISH OPTIONS



Main image shown in **Gray** (standard). **Also available in White** (shown here). Additional frame colors may be available upon request.