

ULTIMATE ERGO MASS SPEC BENCH


WITH DUAL PUMP ENCLOSURE

OMNILAB

SOLUTIONS

OL-MS-ERGO-ULT-D

Overall Bench Dimensions: 96W x 33D x 37H

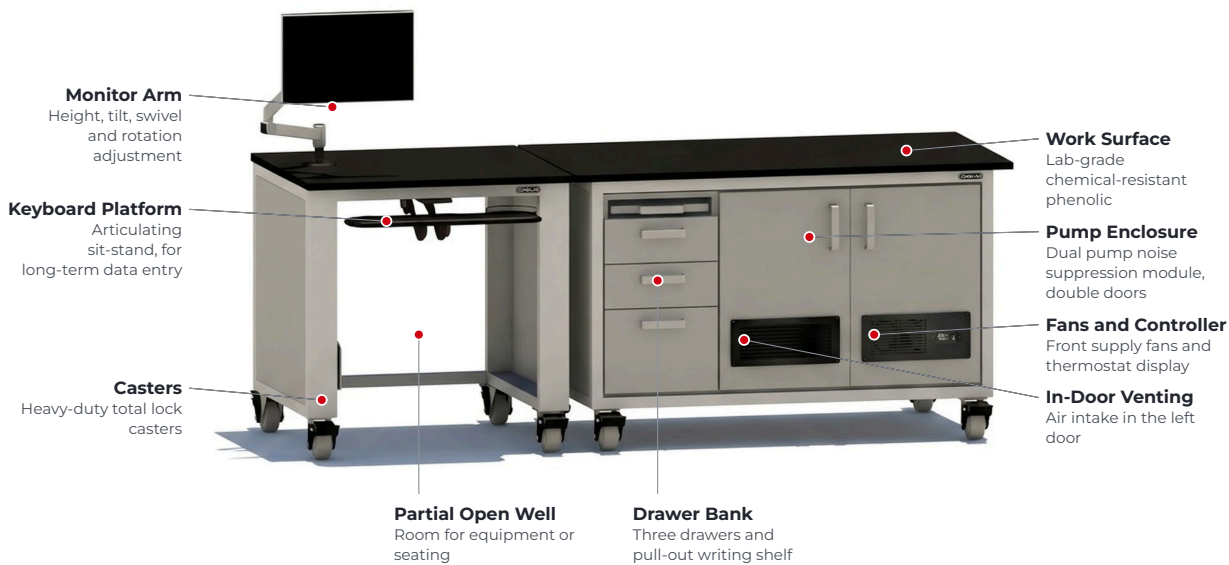

Two-Piece Enclosed Ergo MS Bench (60W MS Bench and 36W open-well section for seating) with lab-grade chemical-resistant phenolic work surfaces and steel frames with lower shelves on heavy-duty total lock casters. Includes a noise suppression module for dual vacuum pumps with acoustic foam paneling, push-pull cooling (440 CFM combined rated exhaust), in-door venting, double doors, five rear hose grommets, and two oil drip pans with shock absorbers; and a drawer bank with three full-extension drawers and a pull-out writing shelf. Also includes a PC accessory kit with a monitor arm, an articulating sit-stand keyboard platform, and a vertical CPU/UPS holder, plus two 120V 8-outlet surge suppressors with 12-foot cords. Overall bench size: 96W x 33D x 37H (Pump Enclosure Interior: 32-1/2W x 25-3/4D x 24H).

SPECIFICATIONS

OVERALL SIZE 96W x 33D x 37H	WORK SURFACE 1" lab-grade chemical-resistant phenolic	FRAME Fully welded steel, epoxy powder coat finish, with wipe-clean laminate shelves
PUMP ENCLOSURE Acoustic foam lined, 440 CFM push-pull cooling, in-door venting, thermostat controller	CASTERS Heavy-duty total lock	STORAGE AND PANELS Wipe-clean laminate in a color that coordinates with the frame

ULTIMATE ERGO MASS SPEC BENCH

WITH DUAL PUMP ENCLOSURE


OMNILAB

SOLUTIONS

COMPONENTS INCLUDED

QTY	MODEL #	DESCRIPTION
1	OL3IT60E	Enclosed Instrument Bench 60W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Inset Steel Frame (48W x 31D x 30-1/4H) with Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity
1	OL3IT36E-FSF-W	Enclosed Instrument Bench 36W x 33D x 37H with Lab-Grade Chemical Resistant Phenolic Work Surface and Flush Steel Frame (35-3/4W x 31D x 30-1/4H) with Partial Open Well for Equipment or Seating and Lower Shelf on Heavy Duty Total Lock Casters, 1,200 lb Uniform Load Capacity
1	OL3NS-D	Bench-Mounted Noise Suppression Module for Dual Vacuum Pumps 36W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with In-Door Venting, Front Supply Fans and Quad Rear Exhaust Fans (440 CFM Combined Rated Exhaust), Five (5) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Double Doors
2	OL3ODP-XL	Oil Drip Pan for Pump with Mounted Shock Absorbers, Large
1	OL3DB16	Drawer Bank 15-3/4W x 22D x 27H with Three (3) Full Extension Lab Drawers and One (1) Pull-Out Utility Shelf
1	OL3MA-C	Ergotron® Work Surface Mounted Ergonomic Monitor Arm, Clamp Mount, Holds Monitors up to 34 Inch and 22 lb, VESA 75 x 75 and 100 x 100 Compatible, Height, Tilt, Swivel and Rotation Adjustment, Integrated Cable Management
1	OL3AKMP	Sit-Stand Articulating Keyboard/Mouse Platform, Stable for Extended Computer Work, 12-1/2 Inch Total Height Adjustment (7 Inch Up, 5-1/2 Inch Down), Dial-In Tilt Adjustment, Fully Retractable Slide Rail Mount
1	OL3CPU	Ergotron® Side-Mounted Vertical CPU/UPS Holder, Adjustable to Fit Units 2 to 8 Inch Wide, 40 lb (18 kg) Capacity, Open Design for Access to Ports and Power Buttons, Steel and Aluminum Construction
2	OL3E	Tripp Lite® Isobar® Surge Suppressor, Eight (8) 120V NEMA 5-15R Outlets in Isolated Filter Banks for EMI/RFI Noise Protection, 3,840 Joule Rating, 12 Amp Resettable Circuit Breaker (1,440 Watt Capacity), UL Listed, Diagnostic LEDs, All-Metal Housing, 12 Ft Cord, Bench Mounted

PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

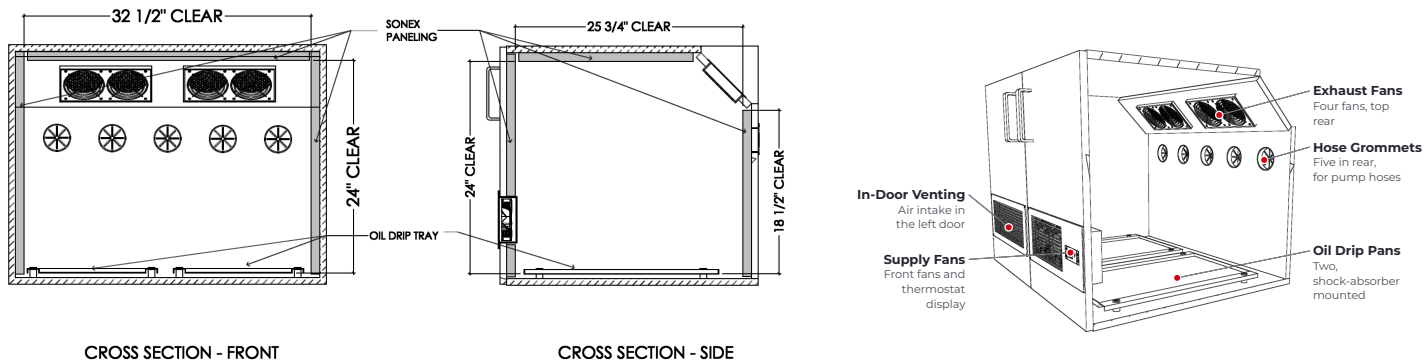

Interior clear dimensions: 32-1/2W x 25-3/4D x 24H (18-1/2H at the rear wall, below the angled exhaust fans). Confirm pumps, hoses, cords and service clearance fit before ordering.

FINISH OPTIONS



Main image shown in **Gray** (standard). **Also available in White** (shown here). Additional frame colors may be available upon request.