


GARAGE-STYLE SINGLE PUMP ENCLOSURE

FOR WHEELED VACUUM PUMPS



OL-NS-G

Overall Dimensions: 22W x 28-1/2D x 27H

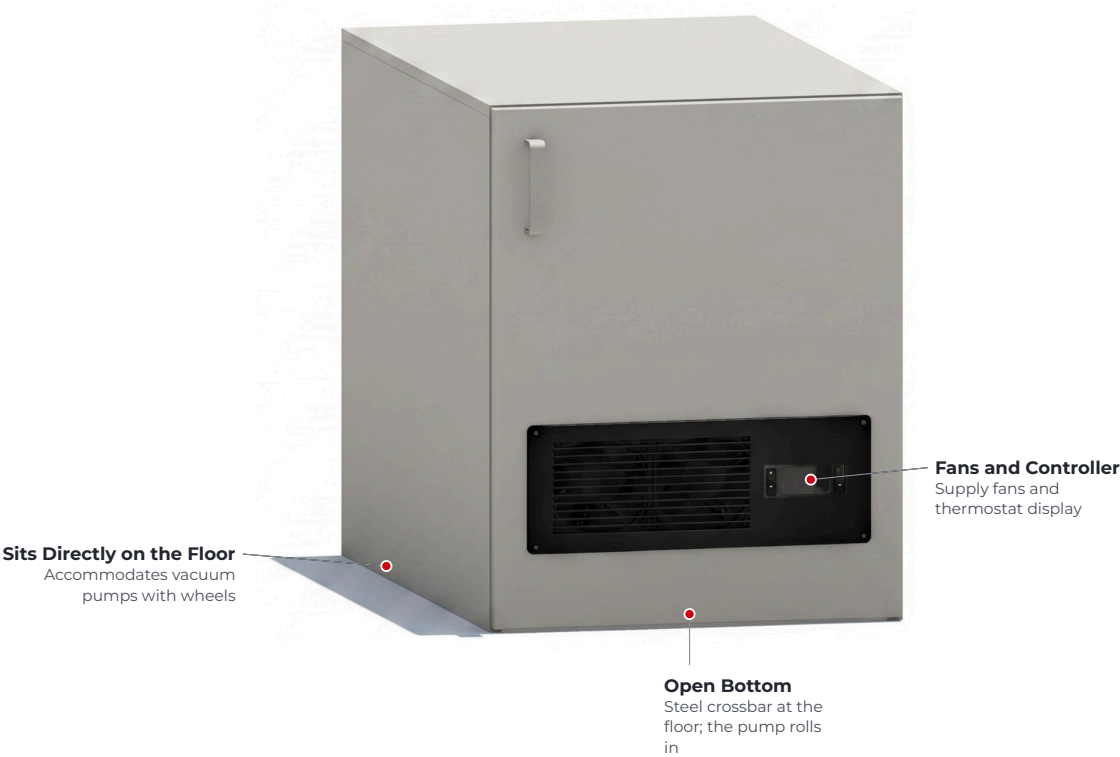

Freestanding Garage-Style Noise Suppression Module for a single wheeled vacuum pump with acoustic foam paneling, push-pull cooling (220 CFM combined rated exhaust), a vented door, three rear hose grommets, and an open bottom with a front steel crossbar that sits directly on the floor so the pump rolls in. Includes a digital display with a programmable thermostat controller, automatic fan speed control, and alarm warnings. Also includes a 6-outlet surge suppressor with a 6-foot cord. Overall size: 22W x 28-1/2D x 27H (Pump Enclosure Interior: 18-1/2W x 25-3/4D x 25H).

SPECIFICATIONS

OVERALL SIZE 22W x 28-1/2D x 27H	PUMP ENCLOSURE Acoustic foam lined, 220 CFM push-pull cooling, thermostat controller	INTERIOR CLEAR 18-1/2W x 25-3/4D x 25H
BASE Open bottom with steel crossbar; no drip pan	CONTROLLER Programmable thermostat, automatic fan speed, alarms	POWER 6-outlet surge suppressor, 6 ft cord

GARAGE-STYLE SINGLE PUMP ENCLOSURE



FOR WHEELED VACUUM PUMPS



COMPONENTS INCLUDED

QTY	MODEL #	DESCRIPTION
1	OL3NSF-G	Freestanding Garage-Style Noise Suppression Module for Single Vacuum Pump 22W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with Dual Front Supply Fans and Dual Rear Exhaust Fans (220 CFM Combined Rated Exhaust), Three (3) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Door Hinged Right, Open Bottom with 1/8 Inch Steel Crossbar at Front, Sits Directly on Floor for Wheeled Pumps
1	OL3E-NS	Tripp Lite® Surge Suppressor, Six (6) 120V NEMA 5-15R Outlets, EMI/RFI Noise Filtering, 2,420 Joule Rating, 15 Amp Circuit Breaker (1,800 Watt Capacity), ETL Listed to UL 1449, Diagnostic LEDs, 6 Ft Cord, Bench Mounted

PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

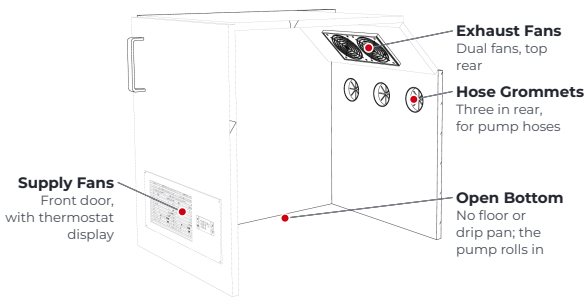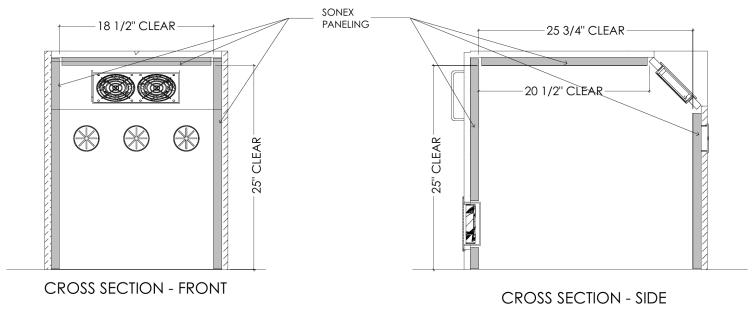

Interior clear dimensions: 18-1/2W x 25-3/4D x 25H (20-1/2D at the top, where the exhaust fans sit at an angle). Open bottom: a wheeled pump rolls straight in. Confirm pump, hoses, cords and service clearance fit before ordering.

FINISH OPTIONS



Main image shown in **Gray** (standard). **Also available in White** (shown here). Additional frame colors may be available upon request.