

SINGLE PUMP ENCLOSURE

WITH FIXED BASE

OL-NS-SF

Overall Dimensions: 22W x 28-1/2D x 31H



Freestanding Noise Suppression Module for a single vacuum pump with acoustic foam paneling, push-pull cooling (220 CFM combined rated exhaust), a vented door, three rear hose grommets, and an oil drip pan with shock absorbers on a fixed base. Includes a digital display with a programmable thermostat controller, automatic fan speed control, and alarm warnings. Also includes a 6-outlet surge suppressor with a 6-foot cord. Overall size: 22W x 28-1/2D x 31H (Pump Enclosure Interior: 18-1/2W x 25-3/4D x 24H).

SPECIFICATIONS

OVERALL SIZE 22W x 28-1/2D x 31H	PUMP ENCLOSURE Acoustic foam lined, 220 CFM push-pull cooling, thermostat controller	INTERIOR CLEAR 18-1/2W x 25-3/4D x 24H
BASE Fixed toe-kick base, 4H	CONTROLLER Programmable thermostat, automatic fan speed, alarms	POWER 6-outlet surge suppressor, 6 ft cord

SINGLE PUMP ENCLOSURE



WITH FIXED BASE



COMPONENTS INCLUDED

QTY	MODEL #	DESCRIPTION
1	OL3NSF-S-R	Freestanding Noise Suppression Module for Single Vacuum Pump 22W x 28-1/2D x 27H with SONEX® Acoustic Foam Paneling on Interior, High-Airflow Push-Pull Cooling System with Dual Front Supply Fans and Dual Rear Exhaust Fans (220 CFM Combined Rated Exhaust), Three (3) Flexible Hose Grommets in Rear, Includes Digital Display, Programmable Thermostat Controller with Automatic Fan Speed Control and Alarm Warnings, Door Hinged Right
1	OL3ODP-XL	Oil Drip Pan for Pump with Mounted Shock Absorbers, Large
1	OL3NSFB-S	Fixed Base for OL3NSF-S 22W x 28-1/2D x 4H
1	OL3E-NS	Tripp Lite® Surge Suppressor, Six (6) 120V NEMA 5-15R Outlets, EMI/RFI Noise Filtering, 2,420 Joule Rating, 15 Amp Circuit Breaker (1,800 Watt Capacity), ETL Listed to UL 1449, Diagnostic LEDs, 6 Ft Cord, Bench Mounted

PUMP ENCLOSURE INTERIOR: CLEAR DIMENSIONS

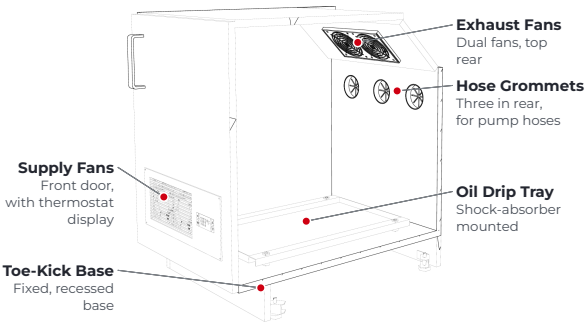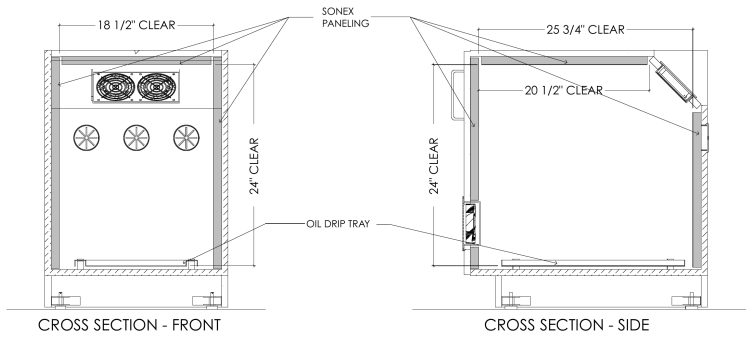

Interior clear dimensions: 18-1/2W x 25-3/4D x 24H (20-1/2D at the top, where the exhaust fans sit at an angle). Confirm pump, hoses, cords and service clearance fit before ordering.

FINISH OPTIONS



Main image shown in **Gray** (standard). **Also available in White** (shown here). Additional frame colors may be available upon request.